
*Escali*

HC255W

**High Capacity Anti-Slip  
Bathroom Scale**

Model #HC255W

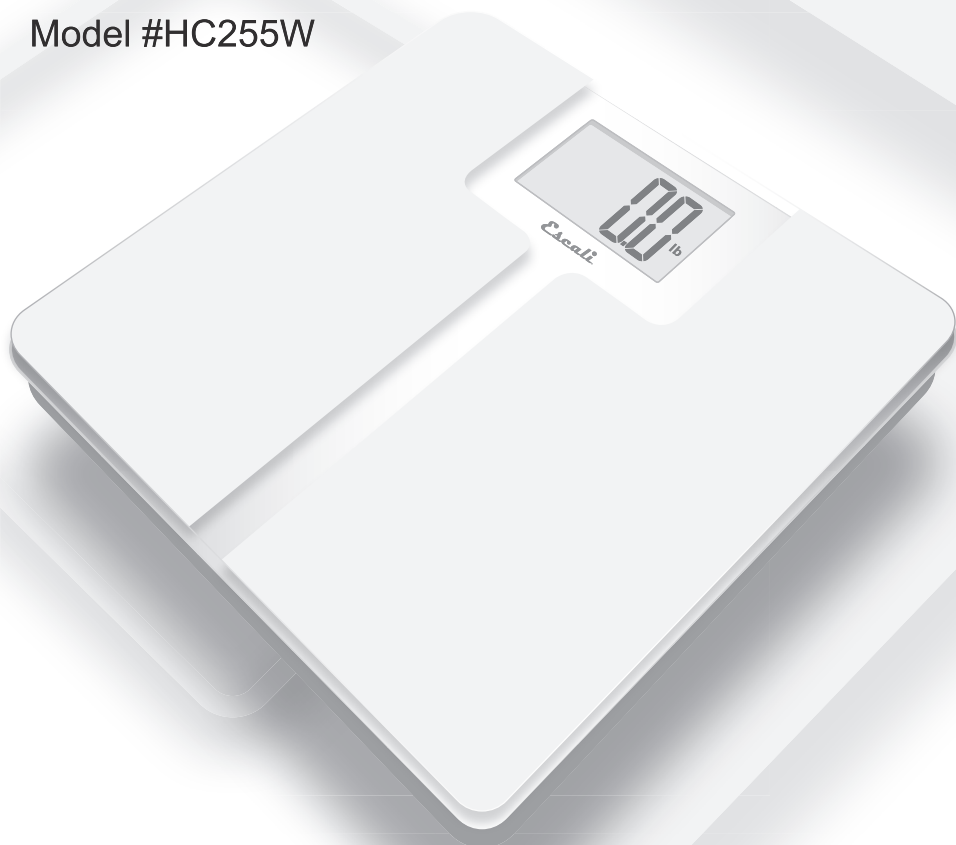
**PRODUCT  
MANUAL**

# High-Capacity Anti-Slip Bathroom Scale

Thank you for purchasing the Escali High-Capacity Anti-Slip Bathroom Scale. Please read the following instructions carefully before using your scale to ensure proper use and care. If you have an issue with your scale, please contact Escali at **800-467-6408** to speak with a product specialist. We would appreciate the opportunity to resolve any issues you may have before returning the product to the retailer.

## Function Diagram

1. Backlit Display
2. Weighing Platform
3. Slip Resistant Etched Glass Surface
4. Unit Button:  
Change between  
lb, kg, and st.
5. Battery Compartment
6. Feet
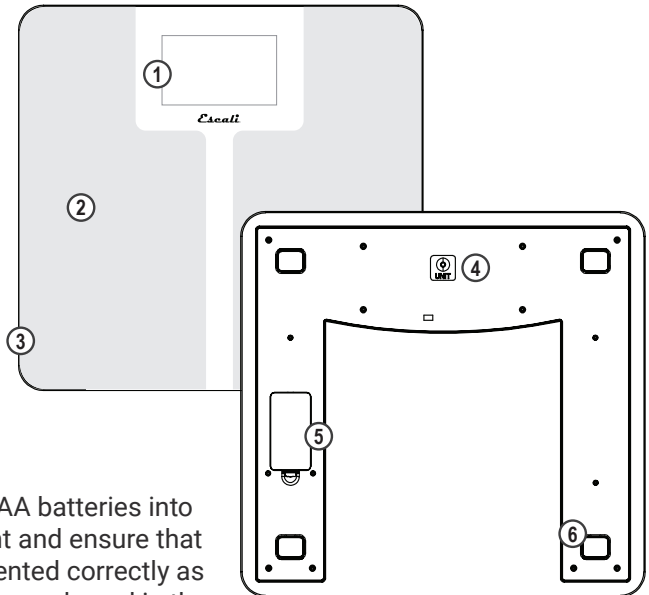
## Operation

1. Place the included 4 x AAA batteries into the battery compartment and ensure that the +/- terminals are oriented correctly as shown on the base of the scale and in the battery compartment.
2. The unit of measure is preset at pounds (lb). This can be changed to kilograms (kg) or stones (st) by pressing the unit button located on the back of the scale.
3. Place the scale on a solid and flat surface.
4. To activate, simply step onto the scale platform. The display will turn on and blink three times when the weight is confirmed. The scale will also produce an audible beep letting you know that the weight has been confirmed.
5. Step off the scale and the weight will remain on screen for 10 seconds before shutting off automatically.

**⚠ NEVER step onto a wet scale or step on the scale platform with wet feet to avoid injury.**

## Scale Recalibration

Your scale is a precision instrument that may need to be recalibrated before use, especially when the scale is picked up or moved. When the scale is moved or jostled and used before being recalibrated, it can result in fluctuations of the weight displayed on screen. To recalibrate the scale, step on the platform to turn the unit on. Step back off of the scale and wait for it to turn back off. This recalibration sequence should remove any inconsistencies and ensure an accurate weight.

**Note:** *The weight displayed on the scale can also fluctuate due to inconsistencies in how you stand on the platform. Stepping on the scale in a consistent manner with your weight evenly distributed will help ensure consistent readings.*

## Safe Use and Care

- Use and stow the scale at room temperature for optimal performance.
- Use on a flat, level, and stable surface to prevent inaccuracy and fluctuation.
- Do Not submerge in liquid. The scale is not waterproof.
- Do Not store in high heat environments or near anything that generates heat. Store in a cool, dry place.
- Do Not use solvents to clean the scale. Use soap and water on a damp sponge or cloth.
- Do Not disassemble or tamper with the scale as it will void the warranty.
- The scale may be affected by electromagnetic disturbances such as radio waves, microwaves, or cell phones. If such disturbances occur, remove the source of the disturbance and restart the weighing process.
- Only install new batteries of the same type. Do Not mix battery types (Alkaline, Carbon-Zinc, rechargeable, etc.) or new and old batteries.
- Recycle or dispose of batteries per state and local guidelines.

 The scale will not weigh properly on carpet.

 The scale is an instrument of precision. Do not submit the scale to shaking, excessive vibration, or other rough treatment as it can permanently damage the scale.

 As the electronic sensors are sensitive, avoid dropping or jarring the scale and store it where it will be protected from impact.

## Warning Messages and Indicators

**EEEE :**

Overload warning – the weight is beyond the scales maximum capacity.

**Lo :**

Low battery warning – the battery power is low and batteries need to be changed.

<b>Product Specifications:</b>	
<b>Max Capacity:</b>	562 lb / 255 kg
<b>Weighing Range:</b>	11 – 562 lb / 5 – 255 kg
<b>Tolerance:</b>	0-50kg ±0.3kg, 50-100kg ±0.4kg, 100-150kg ±0.5kg, 150kg-180kg ±0.6kg, 180-255kg ±1.2kg
<b>Increments:</b>	5-200kg: 0.1 lb, 0.05 kg, 200kg+: 0.1 lb / 0.1 kg
<b>Units:</b>	lb / kg / st
<b>Default Unit:</b>	lb
<b>Response Time:</b>	8 seconds
<b>Weight Confirmation:</b>	Weight blinks three times and audible beep
<b>Battery Required:</b>	Yes (included). Qty. 4 x AAA
<b>Auto Off:</b>	10 seconds
<b>Waterproof:</b>	No
<b>Display:</b>	Backlit LCD (blue)
<b>Display Hold:</b>	10 seconds
<b>Weight:</b>	5.23 lb / 2.37 kg
<b>Dimensions:</b>	12.80 x 12.2 x 1.18" / 325 x 310 x 30mm

## Product Support

24/7 Online support: <https://escali.zendesk.com>

Email support: [feedback@escali.com](mailto:feedback@escali.com)

Phone support: **1-800-467-6408**

## LIFETIME LIMITED WARRANTY

This scale is guaranteed to be free of defects in materials and workmanship.

For more details visit: [www.escali.com/warranty-support](http://www.escali.com/warranty-support)

Escali's liability is limited to two (2) times the cost of the product.

*Escali*

**ESCALI CORPORATION**

3203 Corporate Center Drive, Suite 150

Burnsville, MN 55306

USA

©2021. Escali Corp., Minneapolis, MN, USA. All rights reserved.

Escali, the Escali logo and other Escali marks are owned by Escali Corp. and may be registered.

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© 2021. Escali Corp., Minneapolis, MN, EE. UU.