
*Escali*

**PRODUCT  
MANUAL**

# **Body Analyzing Digital Scale**

Model # BF180

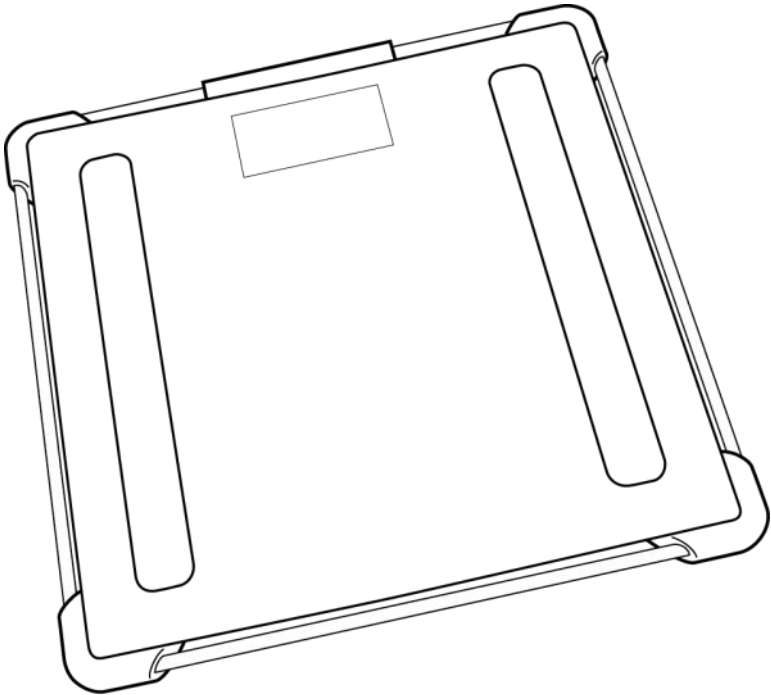
**[www.escali.com](http://www.escali.com)**

Thank you for choosing the Escali BF180 Body Analyzing Digital Scale. This state-of-the-art, precision measuring instrument is only accurate when used on hard, flat and level surfaces. This scale will not function correctly on carpet. The electronic sensors are very sensitive; avoid dropping or jarring the scale, and please keep it where it will be protected from impact.

**Warning!** Do not use this product if you have a pacemaker or other internal electronic medical devices.

## DISPLAY DIAGRAM
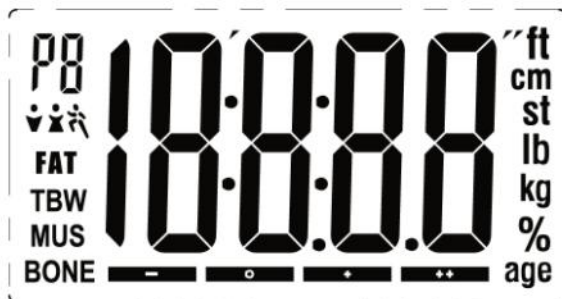
<b>FAT</b>	Body Fat Analysis Result		Male
<b>TBW</b>	Total Body Water Analysis Result		Male Athlete
<b>MUS</b>	Muscle Mass Analysis Result		Woman
<b>BONE</b>	Bone Mass Analysis Result		Woman Athlete
<b>st</b>	Stone		Underfat
<b>lb</b>	Pound		Healthy
<b>kg</b>	Kilogram		Overfat
<b>%</b>	Percentage		Obese
<b>ft</b>	Foot	<b>P8</b>	User Profile Number (P1 to P8)
<b>cm</b>	Centimeter	<b>age</b>	Age

## INITIAL SETUP: CREATING & EDITING USER PROFILES

1. Remove the insulating plastic strip from the battery compartment. Ensure the batteries are properly seated within the +/- terminals.
2. Press the **UNIT** button on the base of the scale to select your desired unit of measure; pounds (lb), kilograms (kg) or stones (st).

- To set up a new user profile (or edit an existing profile) press **SET** and use the **UP** or **DOWN** arrows to select a profile number. Press **SET** to confirm.
- Select your gender using the **UP** or **DOWN** arrows. (Activate Athlete Mode by selecting your gender with the accompanied running figure.) Press **SET** to confirm.
- Select your height using the **UP** or **DOWN** arrows and press **SET** to confirm.
- Select your age using the **UP** or **DOWN** arrows and press **SET** to confirm.
- When the display shows zero weight, step on the scale with clean and dry bare feet. The scale will measure your weight and then calculate your body composition.

**Note:** This body scale can save up to 8 unique user profiles. To program additional profiles, follow steps 3 through 7.

### SHOULD I ACTIVATE ATHLETE MODE?

Activate Athlete Mode if you exercise at least 10 hours per week and have been doing so consistently for the last 6 months, or longer.

### MEASURING BODY WEIGHT & BODY COMPOSITION

- To use the scale, stand still on the center of the platform with clean and dry bare feet touching each of the metal foot placers. After a few seconds the scale will display your weight. 
- Next, the scale will automatically select your user profile and calculate your body composition while a "o" scrolls across the display.

- Finally, the scale will display your current body composition measurements. You may now step off the scale. The scale will cycle through your body fat, body water, muscle mass and bone density measurements three times before itself off.

**Note:** The scale will automatically select your user profile to calculate your body composition by comparing your current body weight measurement to previous body weight measurements from each user profile. If your current body weight changes by more than 7 pounds from your previous measurement, or if your body weight is within 7

pounds of another user profile, then you will need to manually select your user profile number with the corresponding arrow buttons before the scale will display your current body composition.

### HOW DOES THE SCALE MEASURE BODY COMPOSITION?

This scale uses Bioelectrical Impedance Analysis (BIA) technology to measure body composition. BIA technology sends an unperceivable electronic signal through your body to calculate your body fat, body water and muscle mass percentages. Do not use this scale if you have a pacemaker or other internal electronic medical device.

### HOW DO I GET THE MOST ACCURATE RESULTS?

Hydration levels fluctuate throughout the day and are affected by eating, drinking and exercise. As a result, for maximum accuracy and repeatability, it is recommended that this scale is used at the same time of day under the same conditions to establish a personal average. These readings are a close indication, and not an exact value. Please consult your doctor for more information.

### MEASUREMENTS DURING WEIGHT LOSS

During weight loss, body fat percentages may not change as expected since the reading is a percentage of body fat compared to body weight, and both could be changing at the same time. To minimize the loss of lean body tissue, be sure to include exercise in your weight loss plan.

### BODY FAT PERCENTAGE — SUGGESTED RANGES

Age	Men			
	Underfat	Healthy	Overfat	Obese
10-39	< 8%	8-20%	20-25%	25-80%
40-59	< 11%	11-22%	22-28%	28-80%
60-80	< 13%	13-25%	25-30%	30-80%

Age	Women			
	Underfat	Healthy	Overfat	Obese
10-39	< 21%	21-33%	33-39%	39-80%
40-59	< 23%	23-34%	34-40%	40-80%
60-80	< 24%	24-36%	36-42%	42-80%

### MUSCLE MASS PERCENTAGE — SUGGESTED RANGES

Age	6-12	13-20	21-40	41-60	61-80
Men	> 39%	> 37%	> 36%	> 35%	> 33%
Women	> 35%	> 33%	> 32%	> 31%	> 30%

## BODY WATER PERCENTAGE — SUGGESTED RANGES

Age	Men			Women		
	Low	Normal	High	Low	Normal	High
10-15	< 58%	58-72%	> 72%	< 57%	57-67%	> 67%
16-30	< 53%	53-67%	> 67%	< 47%	47-57%	> 57%
31-60	< 47%	47-62%	> 62%	< 42%	42-52%	> 52%
61-80	< 42%	42-56%	> 56%	< 37%	37-47%	> 42%

## BONE MASS (LB/KG) — SUGGESTED RANGES

Body Weight	Men		
	< 132.3 lb	132.3 - 165.3 lb	> 165.3 lb
	< 60 kg	60 - 75 kg	> 75 kg
Est. Bone Mass	5.5 lb	6.5 lb	7 lb
	2.5 kg	2.9 kg	3.2 kg

Body Weight	Women		
	< 99.1 lb	99.1 - 132.3 lb	> 132.3 lb
	< 45 kg	45 - 60 kg	> 60 kg
Est. Bone Mass	4 lb	5 lb	5.5 lb
	1.8 kg	2.2 kg	2.5 kg

## PROPER USE & CARE

- Use and stow the scale at room temperature.
- Do not submit the scale to shaking, excessive vibration, or other rough treatment; use on a flat and stable surface.
- The display may be affected by electromagnetic disturbances such as radios, microwaves, or cell phones. If such disturbances occur, remove the source of the disturbance, and restart the scale.
- Disassembling, or tampering with the scale voids the warranty.
- Clean the scale with a damp cloth; do not submerge in liquid; scale is not waterproof.
- Failure to insert batteries in the correct polarity, as indicated in the battery compartment, may shorten the life of the batteries or cause batteries to leak.
- Do not mix Alkaline, Standard (Carbon-Zinc) or Rechargeable (Nickel Cadmium) or (Nickel Metal Hydride) batteries. Do not mix old and new batteries.
- Batteries should be recycled or disposed of as per state and local guidelines. Do not dispose of batteries in fire.

## TROUBLESHOOTING WARNING MESSAGES

### “- - -” — Overload

*Solution: Stop using the scale for measurement.*

### “Lo”— Low Battery

*Solution: Replace both batteries as measurements may lose accuracy.*

### “Err” — Body Fat Out of Range / BIA Analysis Failure / Sensor Defect

*Solution: Check and confirm your user profile data, stand still on the center of the scale with clean and dry bare feet touching both metal plates and then measure again.*

## PRODUCT SPECIFICATIONS

Model: BF180

Capacity: 400 lb / 180 kg

Increments: 0.2 lb / 0.1 kg

Product Dimensions: 12.875 x 12.875 x 1 in / 30.2 x 30.2 x 2.5 cm

Power Source: Two CR2032 batteries (included)

## ESCALI CORPORATION

3203 Corporate Center Drive, Suite 150

Burnsville, MN 55306

USA

## PRODUCT SUPPORT

24/7 Online support: <https://escali.zendesk.com>

Email support: [feedback@escali.com](mailto:feedback@escali.com)

Phone support: 1-800-467-6408

## LIFETIME LTD. WARRANTY

This scale is guaranteed to be free of defects in materials and workmanship.

For more details visit: [www.escali.com/warranties](http://www.escali.com/warranties)

Escali's liability is limited to two (2) times the cost of the product.

© Escali Corp., Minneapolis, MN, USA. All rights reserved. Escali, the Escali logo and other Escali marks are owned by Escali Corp. and may be registered.

## FCC USER GUIDE INFORMATION

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna;
- Increase the separation between the equipment and the receiver;
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected;
- Consult the dealer or an experienced radio / TV technician for help.

### Radio Transmitters (Part 15)

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

### Caution

The Escali BF180 Body Analyzing Digital Scale offers you a seamless way to manage your health. Please be aware that this device is designed for adults' self-measuring and self-monitoring body fat and hydration levels. Any information provided by this device is in no way meant to treat, cure or prevent any disease or illness from happening. This device should not be used by anyone who is acutely or chronically ill, suffering from a disease or taking medications that affect your water levels. The accuracy of readings for these patients has not been verified. Specific medical advice should be obtained from a physician.

- Before diet and exercise, you must first find a professional medical guide.
- The Escali BF180 is equipped with data transmission function. It may emit electromagnetic energy so as to perform its intended function. Nearby portable and mobile RF communications equipment can affect the performance of the device.
- Portable and mobile RF communications equipment can affect the measuring accuracy of the Escali BF180.
- Kindly note that the use of accessories, transducers or cables other than those specified, with internal components, may result in increased EMISSIONS or decreased IMMUNITY of the device.
- Be aware that misuse of electrical equipments can cause electric shock, burns, fire and other hazards. Please do not use the device in a flammable gas environment.
- Warning that the Escali BF180 should not be used adjacent to or stacked with other equipment.
- Manufacturer will make available on request circuit diagrams, component parts list, etc.
- WARNING: No modifications of this equipment is allowed.
- During using the scale, the patient will contact with the surface and the electrode of the scale. The surface and the electrode of the scale is type BF Applied part.

- The device doesn't need to be calibrated in five years lifetime.
- Do not touch the output of the batteries when using the analyzer.
- The patient is an intended operator. The patient can perform all the operations in the manual, such as measurement, data transmitting, changing batteries.

<b>Product Name</b>	BF180 Body Analyzing Digital Scale
<b>Dimension</b>	Scale: 12.875 x 12.875 x 1 in / 302 x 302 x 25 mm (Approximately)
<b>Net Weight</b>	3.09 lb / 1.40 kg (Approximately)
<b>Display</b>	Digital LED
<b>Measurement Unit</b>	Pound (lb) / Kilogram (kg) / Stone (st)
<b>Measurement Range</b>	11 lb to 397 lb / 5 kg to 180 kg / 0 st: 11 lb to 28 st: 5 lb
<b>Division</b>	0.2 lb / 0.1 kg / 0 st: 1 lb
<b>Accuracy</b>	50kg: ±0.3kg; 100kg: ±0.4kg; 150kg: ±0.5kg; 180kg: ±0.7kg
<b>Working Environment</b>	Temperature: 41°F to 104°F / 5°C to 40°C; Relative Humidity: ≤90% RH; Atmospheric Pressure: 86 kPa to 106 kPa
<b>Storage Environment</b>	Temperature: -4°F to 140°F / -20°C to 60°C; Relative Humidity: ≤90% RH; Atmospheric Pressure: 50 kPa to 106 kPa
<b>Power Source</b>	6V (Two CR2032-size Batteries)
<b>Turn on Method</b>	Instant On Technology
<b>Auto-OFF</b>	The scale will turn off after about 15 seconds if there is no operation.
<b>Accessories</b>	1. Two CR2032-size batteries; 2. User Manual
<b>Mode of Operation</b>	Continuous Operation
<b>Degree of Protection</b>	Type BF applied part
<b>Software Version</b>	A0
<b>Protection against Ingress of Water</b>	IP21

	Symbol for "THE OPERATION GUIDE MUST BE READ"		Symbol for "MANUFACTURER"
	Symbol for "TYPE BF APPLIED PARTS"		Symbol for "DIRECT CURRENT"
<b>SN</b>	Symbol for "SERIAL NUMBER"	 	Symbol for "ENVIRONMENT PROTECTION—Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice."
	Symbol for "MANUFACTURE DATE"		

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