



DETECTO®

OVER **100** YEARS
OF MEDICAL ACCURACY

STAINLESS STEEL DIGITAL SCALE


MODEL D133



USER MANUAL

Congratulations! You have purchased a precision health monitoring product that should provide you with years of accurate use. **This scale has a special finger print /toe print resistant protective coating.** To assure you obtain the best performance please read the following set-up and use instructions carefully before using the scale.

BEFORE USING THIS SCALE—PLEASE READ

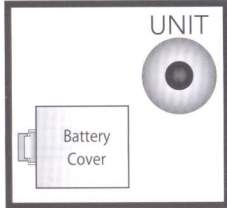


Measurement of Weight

This scale is calibrated to weigh in pounds (lbs), but can also measure in kilograms.

You can convert the weight unit as follows:

1. Find unit switch on the back of the scale.
2. Press the switch to convert scale unit between lb-kg.

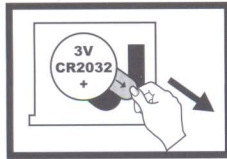

Remove Plastic Tab from Battery

This scale includes one 3V lithium battery (CR2032).

To activate battery: Gently pull out the plastic tab.


OPERATION

1. Place the scale on a flat, dry surface.
2. Gently step on the scale. The scale will automatically turn on.
3. Stand evenly on the scale without moving and wait until your weight, as shown on display, is stable and locked.




AUTOMATIC SWITCH OFF

Auto shut-off occurs when display shows "0.0" for 8 seconds, or shows the same weight reading for 8 seconds.

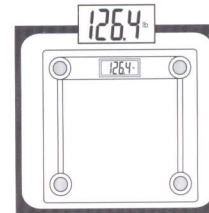


TO REPLACE BATTERY

1. Open battery case cover on the back of scale.
2. Take out the used battery with the aid of a sharp pointed object, as shown in picture to the right.
3. Install new battery by putting one side of the battery beneath the battery baffle and then press down the other side. Be sure (+) side is up.
4. Replace battery cover.

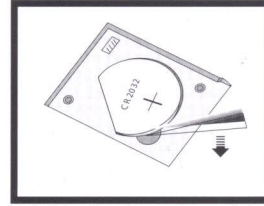


WARNING INDICATIONS

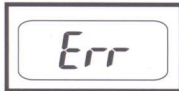

OVERLOAD INDICATION

Scale overload. (Max. weight capacity 400lb/180kg is indicated on the scale). Please step off to avoid damage.

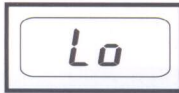

LOW BATTERY INDICATION

Battery power is low. Please install new battery.

ADVICE FOR USE & CARE

1. The stainless steel platform will be slippery when it is wet. Keep dry.
2. Keep still during the course of weighing.
3. Treat the scale with care as it is a precise instrument.
4. Clean the scale with a damp cloth and avoid water getting inside the scale.
DO NOT use chemical/abrasive cleaning agents.
5. Keep the scale in cool place with dry air.
6. If scale fails to turn on, check whether battery is installed properly, or whether battery power is low. If necessary, replace with a new battery.
7. If there is error on display or the scale will not turn off for a long period of time, please take out the battery for about 3 seconds, then install it again to correct the error. If you cannot resolve the problem, please contact Detecto Customer Service for consultation 800-336-5282.
8. Do not use for commercial purposes.

LIMITED LIFETIME WARRANTY

Your DETECTO product is warranted to the original consumer only, for the lifetime of the product, to be free from defects in materials and workmanship under normal use. This product is designed for personal use only. DETECTO, used under license by Lifespan Brands LLC, has the right to nullify the Warranty Agreement should the product be used in a commercial environment. This warranty is void if the product has been tampered with, misused, abused or modified. Our obligation under this Warranty is limited solely to the charge of repairing the product, including the cost of parts, or, at our option, replacing the product with an equivalent model. There are no express or implied warranties except as listed above. DETECTO, USED UNDER LICENSE BY LIFESPAN BRANDS LLC, IS NOT LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES FOR FAILURE OR MALFUNCTION OF THE PRODUCT AND ANY RECOVERY MAY NOT BE GREATER THAN THE PURCHASE PRICE PAID FOR THIS PRODUCT. EXCEPT TO THE EXTENT PROHIBITED BY APPLICABLE LAW, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE ON THIS PRODUCT IS LIMITED TO A PERIOD OF ONE (1) YEAR FROM THE PURCHASE DATE. Some states do not allow limitations on how long an implied warranty lasts, or the exclusion or limitation of incidental or consequential damages, so the above limitations may not apply to you. This warranty gives you specific legal rights and you may have other rights which vary from state to state.

To obtain warranty service, you must return the product freight pre-paid, in either its original packaging or packaging affording an equal degree of protection, to DETECTO, used under license by Lifespan Brands LLC, at 1200 Thorndale Avenue, Elk Grove Village, IL 60007.



For additional information regarding this product or its warranty, please contact us:

DETECTO is a registered trademark of Cardinal Scale Manufacturing Co. and used under license by Lifespan Brands LLC.

Lifespan Brands LLC.
1200 Thorndale Ave., Elk Grove Village, IL 60007.
1-800-336-5282

email: service@detectohome.com

© 2016 Lifespan Brands LLC. All rights reserved.

