
**Ultra Slim  
Health Monitor Scale**

Model # USHM180G-2

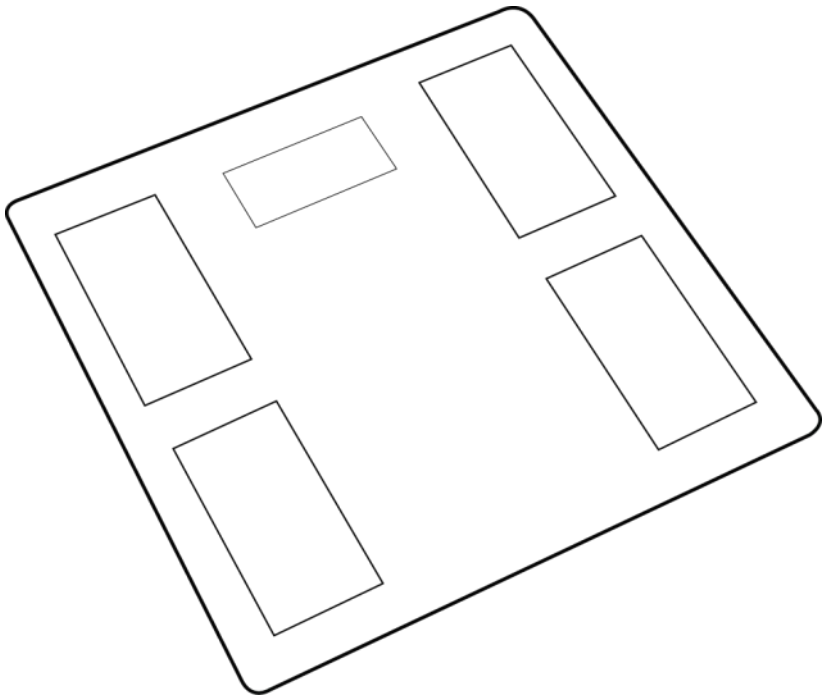
**[www.escali.com](http://www.escali.com)**

Thank you for choosing the Escali USHM180G-2 Goal Tracking Body Composition Scale. This state-of-the-art, precision measuring instrument is only accurate when used on hard, flat and level surfaces. This scale will not function correctly on carpet. The electronic sensors are very sensitive; avoid dropping or jarring the scale, and please keep it where it will be protected from impact.

**Warning!** Do not use this product if you have a pacemaker or other internal electronic medical devices.

## INITIAL SETUP: CREATING & EDITING USER PROFILES

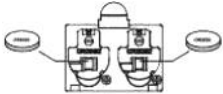1. Remove the insulating plastic strip from the battery compartment. Ensure the batteries are properly seated within the +/- terminals. 
2. Press the 'UNIT' button on the base of the scale to select your desired unit of measure; pounds (lb), kilograms (kg) or stones (st : lb).
3. To set up a new user profile (or edit an existing profile) press **SET** and use the **UP** or **DOWN** arrows to select a user profile number. Press **SET** to confirm.
4. "CLrN" will appear on display. Using the arrow buttons, select "N" to continue, or "Y" to reset your profile. Press **SET** to confirm.
5. Select your gender using the **UP** or **DOWN** arrows. (Activate Athlete Mode by selecting your gender with the accompanied running figure.) Press **SET** to confirm.
6. Select your height using the **UP** or **DOWN** arrows and press **SET** to confirm.
7. Select your age using the **UP** or **DOWN** arrows and press **SET** to confirm.
8. Set your target weight by using the **UP** or **DOWN** arrows and press **SET** to confirm.
9. When the display shows zero weight, step on the scale with clean and dry bare feet. The scale will measure your weight and then calculate your body composition.

**Note:** This body scale can save up to 10 unique user profiles. To program additional profiles, follow steps 3 through 9.

## SHOULD I ACTIVATE ATHLETE MODE?

Activate Athlete Mode if you exercise at least 10 hours per week and have been doing so consistently for the last 6 months, or longer.

## MEASURING BODY WEIGHT & BODY COMPOSITION

1. To use the scale, stand still on the center of the platform with clean and dry bare feet touching each of the metal foot placers. After a few seconds the scale will display your weight.

2. Next, the scale will automatically select your user profile and calculate your body composition while a “o” scrolls across the display.

3. Finally, the scale will display your current body composition measurements. You may now step off the scale. The scale will cycle through your body weight, body fat, body water and muscle mass measurements three times before displaying your user profile information and turning itself off.

**Note:** The scale will automatically select your user profile to calculate your body composition by comparing your current body weight measurement to previous body weight measurements from each user profile. If your current body weight changes by more than 4.5 pounds from your previous measurement, or if your body weight is within 4.5 pounds of another user profile, then you will need to manually select your user profile and confirm your user profile data before the scale will calculate your body composition.

**Note:** If the scale cannot select a user profile automatically it will display for example “P1P2”. You can select the correct profile number by pressing the **UP** or **DOWN** buttons. Press **SET** to confirm, and the scale will proceed to display body composition measurements.

## HOW DOES THE SCALE MEASURE BODY COMPOSITION?

This scale uses Bioelectrical Impedance Analysis (BIA) technology to measure body composition. BIA technology sends an unperceivable electronic signal through your body to calculate your body fat, body water and muscle mass percentages. Do not use this scale if you have a pacemaker or other internal electronic medical device.

## HOW DO I GET THE MOST ACCURATE RESULTS?

Hydration levels fluctuate throughout the day and are affected by eating, drinking and exercise. As a result, for maximum accuracy and repeatability, it is recommended that this scale is used at the same time of day under the same conditions to establish a personal average. These readings are a close indication, and not an exact value. Please consult your doctor for more information.

## REVIEWING USER PROFILE DATA

To view the personal data of your user profile, press the **UP** or **DOWN** arrows with the scale turned off to select a user profile. The scale will automatically cycle through the personal data before displaying zero weight, at which point you could step on the scale to initiate a measurement.

## MEASUREMENTS DURING WEIGHT LOSS

During weight loss, body fat percentages may not change as expected since the reading is a percentage of body fat compared to body weight, and both could be changing at the same time. To minimize the loss of lean body tissue, be sure to include exercise in your weight loss plan.

### BODY FAT PERCENTAGE — SUGGESTED RANGES

Men—

| Age   | Under fat | Healthy | Over fat | Obese |
|-------|-----------|---------|----------|-------|
| ≤ 29  | < 13%     | 13-20%  | 20.1-23% | > 23% |
| 30-39 | < 14%     | 14-21%  | 21.1-24% | > 24% |
| 40-49 | < 16%     | 16-23%  | 23.1-26% | > 26% |
| 50-59 | < 17%     | 17-24%  | 24.1-27% | > 27% |
| 60 ≥  | < 18%     | 18-25%  | 25.1-28% | > 28% |

Women—

| Age   | Under fat | Healthy | Over fat | Obese |
|-------|-----------|---------|----------|-------|
| ≤ 29  | < 19%     | 19-28%  | 28.1-31% | > 31% |
| 30-39 | < 20%     | 20-29%  | 29.1-32% | > 32% |
| 40-49 | < 21%     | 21-30%  | 30.1-33% | > 33% |
| 50-59 | < 22%     | 22-31%  | 31.1-34% | > 34% |
| 60 ≥  | < 23%     | 23-32%  | 32.1-35% | > 35% |

### TOTAL BODY WATER (TBW) PERCENTAGE — SUGGESTED RANGES

|              | Body Fat % Range | Optimal TBW % Range |
|--------------|------------------|---------------------|
|              | Men              | 4 to 14%            |
| 15 to 21%    |                  | 63 to 57%           |
| 22 to 24%    |                  | 57 to 55%           |
| 25% and over |                  | 55 to 37%           |
| Women        | 4 to 20%         | 70 to 58%           |
|              | 21 to 29%        | 58 to 52%           |
|              | 30 to 32%        | 52 to 49%           |
|              | 33% and over     | 49 to 37%           |

### BONE MASS (LB/KG) — SUGGESTED RANGES

| Body Weight    | Men                   |                                |                       | Women                |                               |                       |
|----------------|-----------------------|--------------------------------|-----------------------|----------------------|-------------------------------|-----------------------|
|                | < 132.3 lb<br>< 60 kg | 132.3 - 165.3 lb<br>60 - 75 kg | > 165.3 lb<br>> 75 kg | < 99.1 lb<br>< 45 kg | 99.1 - 132.3 lb<br>45 - 60 kg | > 132.3 lb<br>> 60 kg |
| Est. Bone Mass | 5.5 lb<br>2.5 kg      | 6.5 lb<br>2.9 kg               | 7 lb<br>3.2 kg        | 4 lb<br>1.8 kg       | 5 lb<br>2.2 kg                | 5.5 lb<br>2.5 kg      |

## MUSCLE MASS PERCENTAGE — SUGGESTED RANGES

| Age   | 6-12  | 13-20 | 21-40 | 41-60 | 61-80 |
|-------|-------|-------|-------|-------|-------|
| Men   | > 39% | > 37% | > 36% | > 35% | > 33% |
| Women | > 35% | > 33% | > 32% | > 31% | > 30% |

## PROPER USE & CARE

- Use and stow the scale at room temperature.
- Do not submit the scale to shaking, excessive vibration, or other rough treatment; use on a flat and stable surface.
- The display may be affected by electromagnetic disturbances such as radios, microwaves, or cell phones. If such disturbances occur, remove the source of the disturbance, and restart the scale.
- Disassembling, or tampering with the scale voids the warranty.
- Clean the scale with a damp cloth; do not submerge in liquid; it is not waterproof.
- Failure to insert batteries in the correct polarity, as indicated in the battery compartment, may shorten the life of the batteries or cause batteries to leak.
- Do not mix Alkaline, Standard (Carbon-Zinc) or Rechargeable (Nickel Cadmium) or (Nickel Metal Hydride) batteries. Do not mix old and new batteries.
- Batteries should be recycled or disposed of as per state and local guidelines. Do not dispose of batteries in fire.

## TROUBLESHOOTING WARNING MESSAGES

### “LO”— Low battery indicator

*Solution: Replace the batteries as measurements may lose accuracy.*

### “Err2” — BIA analysis failure / Measurement failure / Overload warning

*Solution: Let the scale turn off. Step back on the center of the platform with clean and dry bare feet with both feet on the metal foot placers.*

### “P1P2” — Unable to automatically select a user profile

*Solution: Use the UP or DOWN arrow buttons to select your user profile.*

### “- . . .” — Unstable surface

*Solution: Ensure the scale is on a hard, flat, level and smooth surface that is free of any vibrations.*

## PRODUCT SPECIFICATIONS

Model: USHM180G-2

Capacity: 400 lb / 180 kg

Increments: 0.2 lb / 0.1 kg

Product Dimensions: 12.25 x 12.25 x 0.83 in (31.1 x 31.1 x 2.1 cm)

Power Source: Two CR2032 batteries (included)

## FCC USER GUIDE INFORMATION

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna;
- Increase the separation between the equipment and the receiver;
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected;
- Consult the dealer or an experienced radio / TV technician for help.

### Radio Transmitters (Part 15)

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

### Caution

The Escali USHM180G Ultra Slim Health Monitor Scale offers you a seamless way to manage your health. Please be aware that this device is designed for adults' self-measuring and self-monitoring body fat and hydration levels. Any information provided by this device is in no way meant to treat, cure or prevent any disease or illness from happening. This device should not be used by anyone who is acutely or chronically ill, suffering from a disease or taking medications that affect your water levels. The accuracy of readings for these patients has not been verified. Specific medical advice should be obtained from a physician.

- Before diet and exercise, you must first find a professional medical guide.
- The Escali USHM180G is equipped with data transmission function. It may emit electromagnetic energy so as to perform its intended function. Nearby portable and mobile RF communications equipment can affect the performance of the device.
- Portable and mobile RF communications equipment can affect the measuring accuracy of the Escali USHM180G.
- Kindly note that the use of accessories, transducers or cables other than those specified, with internal components, may result in increased EMISSIONS or decreased IMMUNITY of the device.
- Be aware that misuse of electrical equipments can cause electric shock, burns, fire and other hazards. Please do not use the device in a flammable gas environment.
- Warning that the Escali USHM180G should not be used adjacent to or stacked with other equipment.
- Manufacturer will make available on request circuit diagrams, component parts list, etc.
- WARNING: No modifications of this equipment is allowed.
- During using the scale, the patient will contact with the surface and the electrode of the

scale. The surface and the electrode of the scale is type BF Applied part.

- The device doesn't need to be calibrated in five years lifetime.
- Do not touch the output of the batteries when using the analyzer.
- The patient is an intended operator. The patient can perform all the operations in the manual, such as measurement, data transmitting, changing batteries.

|  |   |
|--|---|
| <b>Product Name</b>                        | USHM180G-2 Goal Tracking Body Composition Scale   |
| <b>Dimension</b>                           | Scale: 12.25 x 12.25 x 0.83 in / 311 x 311 x 21 mm (Approximately)  |
| <b>Net Weight</b>                          | 3.82 lb / 1.73 kg (Approximately)   |
| <b>Display</b>                             | LCD   |
| <b>Measurement Unit</b>                    | Pound (lb) / Kilogram (kg) / Stone (st)   |
| <b>Measurement Range</b>                   | 11 lb to 397 lb / 5 kg to 180 kg / 0 st: 11 lb to 28 st: 5 lb   |
| <b>Division</b>                            | 0.2 lb / 0.1 kg / 0 st: 0.2 lb, 1 lb (≥20 st)   |
| <b>Accuracy</b>                            | 50kg: ±0.3kg; 100kg: ±0.4kg; 150kg: ±0.5kg; 180kg: ±0.7kg   |
| <b>Working Environment</b>                 | Temperature: 50°F to 104°F / 10°C to 40°C; Relative Humidity: 15% RH to 90% RH; Atmospheric Pressure: 80 kPa to 105 kPa |
| <b>Storage Environment</b>                 | Temperature: -4°F to 131°F / -20°C to 55°C; Relative Humidity: ≤90% RH  |
| <b>Power Source</b>                        | 6V (Two CR2032-size Batteries)  |
| <b>Turn on Method</b>                      | Instant On Technology   |
| <b>Auto-OFF</b>                            | The scale will turn off after about 15 seconds if there is no operation.  |
| <b>Accessories</b>                         | 1. Two CR2032-size batteries; 2. User Manual  |
| <b>Mode of Operation</b>                   | Continuous Operation  |
| <b>Degree of Protection</b>                | Type BF applied part  |
| <b>Software Version</b>                    | A0  |
| <b>Protection against Ingress of Water</b> | IP21  |

|  |   |   |  |
|--|---|---|--|
|  | Symbol for "THE OPERATION GUIDE MUST BE READ" |   |  |
|  | Symbol for "TYPE BF APPLIED PARTS"            |  | Symbol for "DIRECT CURRENT"  |
|  | Symbol for "SERIAL NUMBER"                    |  | Symbol for "ENVIRONMENT PROTECTION—Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice." |
|  | Symbol for "MANUFACTURE DATE"                 |   |  |

## EMC GUIDANCE

Table 1 —

| <b>Guidance and Manufacturer's Declaration — Electromagnetic Emissions</b> |                   |
|--|-------------------|
| <b>Emissions test</b>  | <b>Compliance</b> |
| RF emissions<br>CISPR 11   | Group 1           |
| RF emissions<br>CISPR 11   | Class B           |
| Harmonic emissions<br>IEC 61000-3-2  | Not application   |
| Voltage fluctuations/ flicker emissions<br>IEC 61000-3-3                   | Not application   |

Table 2 —

| <b>Guidance and Manufacturer's Declaration — Electromagnetic Immunity</b>                                    |  |  |
|--|--|--|
| <b>Immunity Test</b>   | <b>IEC 60601-1-2<br/>Test level</b>              | <b>Compliance level</b>                          |
| Electrostatic discharge (ESD)<br>IEC 61000-4-2   | ±8 kV contact<br>±2 kV, ±4 kV, ±8 kV, ±15 kV air | ±8 kV contact<br>±2 kV, ±4 kV, ±8 kV, ±15 kV air |
| Electrical fast transient/burst<br>IEC 61000-4-4   | Not application                                  | Not application                                  |
| Surge<br>IEC 61000-4-5   | Not application                                  | Not application                                  |
| Voltage dips, short interruptions<br>and voltage variations on power<br>supply input lines<br>IEC 61000-4-11 | Not application                                  | Not application                                  |
| Power frequency magnetic field<br>IEC 61000-4-8  | 30 A/m<br>50Hz/60Hz                              | 30 A/m<br>50Hz/60Hz                              |
| Conducted RF<br>IEC61000-4-6   | Not application                                  | Not application                                  |
| Radiated RF<br>IEC61000-4-3  | 10 V/m<br>80 MHz – 2,7 GHz<br>80 % AM at 1 kHz   | 10 V/m<br>80 MHz – 2,7 GHz<br>80 % AM at 1 kHz   |
| NOTE U <sub>T</sub> is the a.c. mains voltage prior to application of the test level.                        |  |  |

Table 3 —

| Guidance and Manufacturer's Declaration - Electromagnetic Immunity   |                      |             |   |                                 |                |              |                           |
|--|----------------------|-------------|---|---------------------------------|----------------|--------------|---------------------------|
| Radiated RF<br>IEC61000-4-3<br>(Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment) | Test Frequency (MHz) | Band (MHz)  | Service   | Modulation                      | Modulation (W) | Distance (m) | IMMUNITY TEST LEVEL (V/m) |
|  | 385                  | 380 – 390   | TETRA 400   | Pulse modulation 18 Hz          | 1,8            | 0.3          | 27                        |
|  | 450                  | 430-470     | GMRS 460, FRS 460   | FM ± 5 kHz deviation 1 kHz sine | 2              | 0.3          | 28                        |
|  | 710                  | 704 – 787   | LTE Band 13, 17   | Pulse modulation 217 Hz         | 0,2            | 0.3          | 9                         |
|  | 745                  |             |   |                                 |                |              |                           |
|  | 780                  |             |   |                                 |                |              |                           |
|  | 810                  | 800 – 960   | GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band            | Pulse modulation 18 Hz          | 2              | 0.3          | 28                        |
|  | 870                  |             |   |                                 |                |              |                           |
|  | 930                  |             |   |                                 |                |              |                           |
|  | 1720                 | 1700 – 1990 | GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS | Pulse modulation 217 Hz         | 2              | 0.3          | 28                        |
|  | 1845                 |             |   |                                 |                |              |                           |
|  | 1970                 |             |   |                                 |                |              |                           |
|  | 2450                 | 2400 – 2570 | Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7            | Pulse modulation 217 Hz         | 2              | 0.3          | 28                        |
|  | 5240                 | 5100 – 5800 | WLAN 802.11 a/n   | Pulse modulation 217 Hz         | 0,2            | 0.3          | 9                         |
|  | 5240                 |             |   |                                 |                |              |                           |
| 5785   |                      |             |   |                                 |                |              |                           |

## **ESCALI CORPORATION**

3203 Corporate Center Drive, Suite 150  
Burnsville, MN 55306  
USA

### **PRODUCT SUPPORT**

24/7 Online support: <https://escali.zendesk.com>

Email support: [feedback@escali.com](mailto:feedback@escali.com)

Phone support: 1-800-467-6408

### **LIFETIME LTD WARRANTY**

This scale is guaranteed to be free of defects in materials and workmanship. For more details visit: [www.escali.com/warranties](http://www.escali.com/warranties)

Escali's liability is limited to two (2) times the cost of the product.

© Escali Corp., Minneapolis, MN, USA. All rights reserved. Escali, the Escali logo and other Escali marks are owned by Escali Corp. and may be registered.

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### **GARANTIE À VIE LTD.**