T180

**COMPACT**  
BODY SCALE

Model #T180

22i


**PRODUCT  
MANUAL**

# Escali Compact Body Scale

Thank you for purchasing the Escali Compact Body Scale. Please read the following instructions carefully before using your scale to ensure proper use and care. If you have an issue with your scale, please contact Escali at 1-800-467-6408 to speak with a product specialist. We would appreciate the opportunity to resolve any issues you may have before returning the product to the retailer.

## Function Diagram

1. Display
2. Weighing Platform
3. Unit Button
4. Battery Compartment
5. Feet
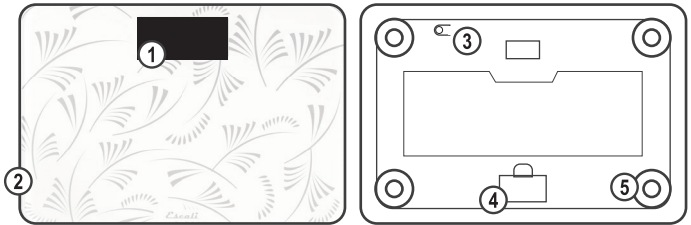
## Operation

Before using your scale for the first time, remove the battery insulating tab which is placed between one of the batteries and the battery terminal in the battery compartment to prevent battery drain during shipping.

1. Ensure that the included 2 x CR2032 batteries are inserted into the battery compartment and oriented correctly with the positive terminal facing up.
2. The unit of measure is preset at pounds (lb). This can be changed to kilograms (kg) or stones (st) by pressing the unit button located on the back of the scale.
3. Place the scale on a solid and flat surface.
4. To activate, simply step onto the scale platform. The display will turn on and blink three times when the weight is confirmed.
5. Step off the scale and the weight will remain on screen for 10 seconds before shutting off automatically.

**⚠ NEVER step onto a wet scale or step on the scale platform with wet feet to avoid injury.**

## Scale Recalibration

Your scale is a precision instrument that may need to be recalibrated before use, especially when the scale is picked up or moved. When the scale is moved or jostled and used before being recalibrated, it can result in fluctuations of the weight displayed on screen. To recalibrate the scale, step on the platform to turn the unit on. Step back off of the scale and wait for it to turn back off. This recalibration sequence should remove any inconsistencies and ensure an accurate weight.

**Note:** The weight displayed on the scale can also fluctuate due to inconsistencies in how you stand on the platform. Stepping on the scale in a consistent manner with your weight evenly distributed will help ensure consistent readings.

## Safe Use and Care

- Use and stow the scale at room temperature for optimal performance.
- Use on a flat, level, and stable surface to prevent inaccuracy and fluctuation.
- Do Not submerge in liquid. The scale is not waterproof.
- Do Not store in high heat environments or near anything that generates heat. Store in a cool, dry place.
- Do Not use solvents to clean the scale or cover. Use soap and water on a damp sponge or cloth.
- Do Not disassemble or tamper with the scale as it will void the warranty.
- The scale may be affected by electromagnetic disturbances such as radio waves, microwaves, or cell phones. If such disturbances occur, remove the source of the disturbance and restart the weighing process.
- Only install new batteries of the same type. Do Not mix battery types (Alkaline, Carbon-Zinc, rechargeable, etc.) or new and old batteries.
- Recycle or dispose of batteries per state and local guidelines.

 The scale will not weigh properly on carpet.

 The scale is an instrument of precision. Do not submit the scale to shaking, excessive vibration, or other rough treatment as it can permanently damage the scale.

 As the electronic sensors are sensitive, avoid dropping or jarring the scale and store it where it will be protected from impact.

## Warning Messages and Indicators

### **ERR :**

Overload warning – the weight is beyond the scales maximum capacity.

### **Lo :**

Low battery warning – the battery power is low and batteries need to be changed.

<b>Product Specifications:</b>	
<b>Max Capacity:</b>	400 lb / 180 kg
<b>Weighing Range:</b>	11 – 400 lb / 5 – 180 kg
<b>Tolerance:</b>	0-50kg ±0.3kg, 50-100kg ±0.4kg, 100-150kg ±0.5kg, 150kg-180kg ±0.7kg
<b>Increments:</b>	0.2 lb / 0.1 kg
<b>Units:</b>	lb / kg / st
<b>Default Unit:</b>	lb
<b>Response Time:</b>	3-5 seconds
<b>Weight Confirmation:</b>	Weight blinks three times
<b>Battery Required:</b>	Yes (included)
<b>Battery Type:</b>	2 x CR2032
<b>Auto Off:</b>	10 seconds
<b>Waterproof:</b>	No
<b>Display:</b>	Backlit LCD
<b>Display Hold:</b>	10 seconds
<b>Weight:</b>	2.68 lb / 1.2 kg
<b>Dimensions:</b>	11.8 x 8.7 x 0.8" / 300 x 220 x 20mm

## Product Support

24/7 Online support: <https://escali.zendesk.com>

Email support: [feedback@escali.com](mailto:feedback@escali.com)

Phone support: **1-800-467-6408**

## LIFETIME LIMITED WARRANTY

This scale is guaranteed to be free of defects in materials and workmanship.

For more details visit: [www.escali.com/warranty-support](http://www.escali.com/warranty-support)

Escali's liability is limited to two (2) times the cost of the product.

***Escali***

**ESCALI CORPORATION**

3203 Corporate Center Drive, Suite 150

Burnsville, MN 55306

USA

©2022. Escali Corp., Minneapolis, MN, USA. All rights reserved.

Escali, the Escali logo and other Escali marks are owned by Escali Corp. and may be registered.

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© 2022. Escali Corp., Minneapolis, MN, EE. UU.